nature neuroscience

Building three-dimensional human brain organoids

Sergiu P. Pasca

The organogenesis of the human central nervous system is an intricately orchestrated series of events that occurs over several months and ultimately gives rise to the circuits underlying cognition and behavior. There is a pressing need for developing reliable, realistic, and personalized in vitro models of the human brain to advance our understanding of neural development, evolution, and disease. Pluripotent stem cells have the remarkable ability to differentiate in vitro into any of the germ layers and, with the advent of three-dimensional (3D) cell

culture methods, to self-organize into brain spheroids or organoids. These organoid cultures can be derived from any individual, can be guided to resemble specific brain regions, and can be employed to model complex cell-cell interactions in assembloids and to build human circuits. This emerging technology, in combination with bioengineering and other state-of-the-art methods for probing and manipulating neural tissue, has the potential to bring insights into human brain organogenesis and the pathogenesis of neurological and psychiatric disorders.



Scientists Helping Scientists[™] | WWW.STEMCELL.COM

Brain organogenesis in vitro and in vivo

Methods for generating neural cells in vitro aim to recapitulate key stages of in vivo brain organogenesis. Folding of the ectoderm-derived neural plate gives rise to the neural tube, which becomes enlarged on the anterior side to form the forebrain in the central nervous system (CNS). Corticogenesis involves the sequential generation and positioning of layerspecific glutamatergic neurons from progenitors that line the ventricles in the dorsal forebrain, the migration of GABAergic interneurons that are born in the ventral forebrain, and waves of gliogenesis to form astrocytes and oligodendrocytes, which continue postnatally.

Pluripotent stem cells, derived from the inner mass of the blastocyst (embryonic stem cells) or from reprogrammed somatic cells (induced pluripotent stem cells), can be differentiated into neural cells in bi-dimensional (2D) cultures - where early on neuroepithelial cells position themselves into structures called rosettes — or in self-organizing threedimensional (3D) brain spheroids or organoids¹. Intermediate ('2.5D') cultures can be obtained when neural cells differentiated in 2D are lifted and cultured in 3D conditions to form cellular aggregates or when differentiated 3D aggregates derived from pluripotent stem cells are subsequently plated for

Brain organoids can be generated from aggregates of pluripotent stem cells through undirected differentiation methods that lack inductive signals^{2,3}, or by patterning through directed differentiation methods to resemble specific brain regions⁴⁻⁸ (e.g., forebrain, midbrain, retina).


Brain assembloids

To model interactions between brain regions, organoids can be patterned to resemble specific regions of the nervous system and then can be fused to generate brain assembloids¹. Another way to generate assembloids is by spatio-temporally controlling patterning within one 3D aggregate, by embedding organizerlike structures (i.e., cells, coated beads) that release or block developmental signals. A third method involves combining other single cells into brain organoids: for instance, by embedding yolk-sac-derived microglia to study neuroimmune interactions, by embedding mesoderm-derived vascular cells to study the blood-brain barrier, or by embedding tumor cells to study brain metastasis. These 3D cultures can be probed using genetic, anatomical and functional read-outs.


Applications of brain organoids

Brain organoids and assembloids can be used to ask questions about evolutionary innovation in human and non-human primates and to understand the developmental program and maturation of the nervous system (e.g., the programs underlying astrocyte transition from a fetal to a postnatal state). Brain organoids derived from patients or that have been genetically engineered to carry genetic variants associated with disease (i.e., isogenic lines) can be used to investigate disease pathogenesis in the nervous system. Lastly, as these 3D cultures become more scalable and assays probing 3D tissue improve, drug and CRISPR-Cas9-based screens can be used to identify therapeutic targets.

Brain organogenesis in vitro and in vivo

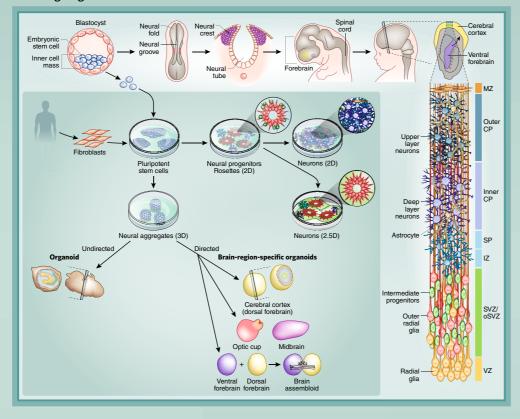

Brain assembloids

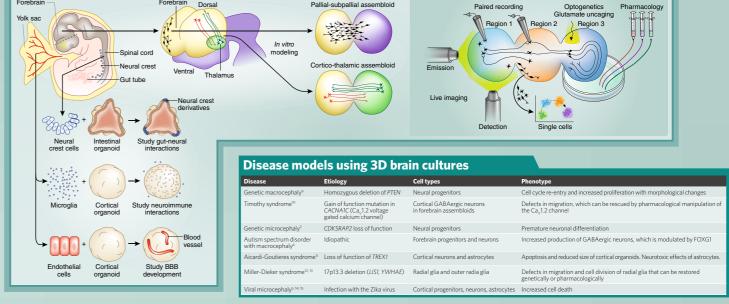

Challenges and future directions

- Improving reliability, anatomical accuracy (cortex expansion, folding, white matter), predictability, and scalability of brain organoids.
- Developing methods for transplantation of organoids into rodents or other species to obtain circuit-wide integration and oscillatory activity, to study sensory input, and to develop more accurate models of


Culturing cells in 3D

psychiatric disease.

Embedding

- Modeling advanced stages, including postnatal, of human brain development and incorporating missing cell types, e.g. glia, endothelial cells and pericytes, to study neuroimmune and neurovascular interactions or cancer cells (e.g., oncogenesis and metastasis). Developing reliable models of environmental,
- monogenic, and polygenic causes of CNS disease to explore questions about convergent and divergent pathogenesis in psychiatric disorders (e.g., autism spectrum disorders, schizophrenia).
- Building large-scale platforms for drug discovery and genetic screens.

Applications of brain organoids

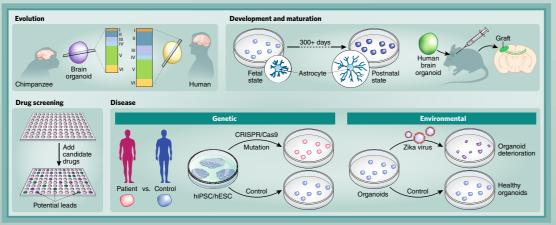


Maintaining 3D cell cultures







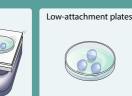

Detachment from plates

Methods for culturing cell aggregation by centrifugation in U-

extracellular matrices and detachment of intact pluripotent colonies that are then moved to ultra-low-attachment dishes. These 3D cultures can be subsequently maintained in low-or high-oxygen conditions, shaker or spun in a bioreactor, or maintained in low-attachment plates without shaking.

STEMCELL Technologies

STEMdiff™ Cerebral Organoid Kit is a 3D in vitro culture system designed to generate cerebral organoids from human pluripotent stem cells. The resulting cerebral organoids have a cellular composition and structural organization representative of the developing human brain.

See the data, at www.stemcell.com/COKit

- STEMdiff™ Cerebral Organoid Kit (#08570) Recapitulates the developmental processes and making sure your research works.
- organization of the developing human brain Optimized based on the formulation published by MA Lancaster and JA Knoblich¹ www.stemcell.com/CerebralOrganoids Rigorous raw material screening and quality control testing ensure reproducibility and
- minimal lot-to-lot variability Simple, serum-free and easy to use

At STEMCELL, science is our foundation. We

BBB: blood-brain barrier; CP: cortical plate; ENS: enteric nervous system; IZ: intermediate zone; MGE: medial nic eminence; MZ: marginal zone; SP: subplate

Aggregation in

U- or V-bottomed wells

Sergiu P. Paşca is at the Department of Psychiatry and Behavioral Sciences, Stanford University, Stanford, CA, USA

© 2018 Nature America, Inc., part of Springer Nature www.nature.com/articles/10.1038/s41593-018-0107-3