SCIENTIFIC REPORTS

Received: 22 June 2018 Accepted: 31 January 2019 Published online: 28 February 2019

OPEN GPDPLQ₁₂₃₇—A Type II Collagen **Neo-Epitope Biomarker of Osteoclast- and Inflammation-Derived Cartilage Degradation** in vitro

Henrik Löfvall^{1,2}, Anna Katri^{1,3}, Aneta Dabrowska¹, Morten A. Karsdal¹, Yunyun Luo¹, Yi He¹, Tina Manon-Jensen¹, Morten H. Dziegiel⁴, Anne-Christine Bay-Jensen¹, Christian S. Thudium¹ & Kim Henriksen¹

C-telopeptide of type II collagen (CTX-II) has been shown to be a highly relevant biomarker of cartilage degradation in human rheumatic diseases, if measured in synovial fluid or urine. However, serum or plasma CTX-II have not been demonstrated to have any clinical utility to date. Here, we describe the GPDPLQ₁₂₃₇ ELISA which targets the EKGPDPLQ↓ neo-epitope, an elongated version of the CTX-II neoepitope (EKGPDP \downarrow), speculated to be a blood-precursor of CTX-II generated by the cysteine protease cathepsin K. Human osteoclast cartilage resorption cultures as well as oncostatin M and tumour necrosis factor α -stimulated bovine cartilage explant cultures were used to validate GPDPLQ₁₂₃₇ biologically by treating the cultures with the cysteine protease inhibitor E-64 and/or the matrix metalloproteinase (MMP) inhibitor GM6001 to assess the potential contributions of these two protease classes to GPDPLQ₁₂₃₇ release. Cartilage resorption-derived GPDPLQ₁₂₃₇ release was inhibited by E-64 (72.1% inhibition), GM6001 (75.5%), and E-64/GM6001 (91.5%), whereas CTX-II release was inhibited by GM6001 (87.0%) but not by E-64 (5.5%). Cartilage explant GPDPLQ₁₂₃₇ and CTX-II release were both fully inhibited by GM6001 but were not inhibited by E-64. No clinically relevant GPDPLO₁₂₃₇ reactivity was identified in human serum, plasma, or urine from healthy donors or arthritis patients. In conclusion, the GPDPLQ₁₂₃₇ biomarker is released during osteoclast-derived cysteine protease- and MMP-mediated cartilage degradation in vitro, whereas CTX-II release is mediated by MMPs and not by cysteine proteases, as well as from MMP-mediated cartilage degradation under a pro-inflammatory stimulus. These findings suggest that GPDPLO₁₂₃₇ may be relevant in diseases with pathological osteoclast activity and cartilage degradation. Further studies are required to validate the neo-epitope in human samples.

C-telopeptide of type II collagen (CTX-II) is a biomarker of type II collagen degradation, indicative of cartilage degradation. Degradation of joint extracellular matrices (ECMs), including cartilage, is a key feature of multiple types of arthritis, including osteoarthritis (OA) and rheumatoid arthritis (RA). Urinary and synovial fluid (SF) CTX-II, unlike serum CTX-II, is a highly relevant biomarker for OA¹⁻³ with some relevance for RA⁴. OA incidence and progression^{5,6}, severity^{6,7-8}, bone marrow lesions⁹, osteophytes¹⁰ and pain¹¹ have been associated with CTX-II, although pain associations have varied¹¹. In RA, CTX-II has been shown to predict disease progression¹²⁻¹⁵, treatment efficacy¹³, bone mineral density (BMD) reduction¹⁶ and synovitis¹⁶. However, no studies have been published to date that demonstrate any clinical utility of serum or plasma CTX-II, although it has been demonstrated to be very useful in animal models with experimentally induced joint destruction¹⁷⁻²¹. Additionally,

¹Nordic Bioscience, Herley, Denmark. ²Division of Molecular Medicine and Gene Therapy, Lund Strategic Center for Stem Cell Biology, Lund, Sweden. ³Department of Drug Design and Pharmacology, Copenhagen University, Copenhagen, Denmark. ⁴Department of Clinical Immunology, Rigshospitalet, Copenhagen University Hospital, Copenhagen, Denmark. Correspondence and requests for materials should be addressed to K.H. (email: kh@ nordicbio.com)

Peptide identifier	Amino acid sequence
Biotinylated peptide	biotin-EKGPDPLQ
Standard peptide	EKGPDPLQ
Elongated +1 peptide	EKGPDPLQY
Elongated +2 peptide	EKGPDPLQYM
Truncated -1 peptide	EKGPDPL
Truncated -2 peptide	EKGPDP
Immunogenic peptide	KLH-CEKGPDPLQ

Table 1. Synthetic peptides used for development and validation of the GPDPLQ1237 competitive ELISA.Abbreviations: KLH, keyhole limpet haemocyanin.

.....

SF and urine sampling have their own drawbacks, particularly in terms of practicality. SF collection is more complicated and inconvenient than collecting either serum, plasma, or urine and is also a more complex biological matrix to work with. Furthermore, both urine and SF analyte concentrations require volume adjustment due to e.g. fluid intake, arthroscopic lavage or joint effusion^{3,22}. Urine CTX-II is also known to have diurnal variation²³, which necessitates a strict urine sampling procedure. As blood is more frequently sampled than urine and SF, a blood-based CTX-II-like biomarker would be more accessible, practical, and reliable for assessing cartilage degradation, if clinical utility could be demonstrated.

CTX-II assays detect C-telopeptide fragments of type II collagen, the major ECM component of articular cartilage, with a neo-epitope²⁴ formed by proteolytic cleavage (EKGPDP \downarrow)²⁵ during degradation of cartilage by MMPs²⁶⁻²⁸. The neo-epitope has been immunolocalised to areas of cartilage damage, around chondrocytes, areas of vascularization close to the subchondral bone, and the bone–cartilage interface²⁹—suggesting possible contribution of osteoclasts to CTX-II. Osteoclasts play important roles in diseases with progressive joint destruction, particularly for bone erosion in diarthrodial joints in RA³⁰⁻³⁴. However, their role in cartilage alterations—such as those seen in OA—is poorly understood, albeit strongly implicated³⁵. An elongated CTX-II neo-epitope (EKGPDPLQ \downarrow) has been previously described in patent literature^{26–28}, where it has been suggested to be generated by cathepsin K-mediated cleavage of calcified cartilage, to be osteoclast-specific, and to be present in blood while being absent or only present in trace amounts in urine and SF. In the patents, proteases in the kidneys and/ or liver were described to potentially remove the C-terminal -LQ of the EKGPDPLQ neo-epitope to generate the urinary forms of the C-telopeptide—EKGPDP or smaller^{26–28}. Due to potential differences in generation and/ or processing from the EKGPDP neo-epitope, the EKGPDPLQ neo-epitope may be an interesting blood-based biomarker for studies of degenerative joint diseases.

To date no peer-reviewed clinical, *in vivo*, or *in vitro* studies using assays detecting the EKGPDPLQ neo-epitope have been published. In this study we investigated the release of this neo-epitope from non-calcified articular cartilage, in *in vitro* models of cartilage degradation derived from osteoclasts³⁶ or inflammation²⁴, using the novel competitive GPDPLQ₁₂₃₇ enzyme-linked immunosorbent assay (ELISA). The validated assay was then used to test for GPDPLQ₁₂₃₇ reactivity in human serum, plasma, and urine.

Results

Specificity of the GPDPLO₁₂₃₇ **ELISA.** The target sequence of the GPDPLQ₁₂₃₇ ELISA, $_{1230}$ EKGPDPLQ₁₂₃₇, was analysed for homology to other human and animal proteins using the NPS@: Network Protein Sequence Analysis PattInProt search against the UniProtKB/Swiss-prot database. The target sequence was found to be unique to the type II collagen alpha 1 chain, and was fully conserved in human, cow, and rat. When allowing for one mismatch in the sequence, it was also conserved in mouse but not in any other animal proteins.

The specificity of the competitive GPDPLQ₁₂₃₇ ELISA was evaluated by analysing the reactivity towards the standard peptide as well as two elongated and two truncated standard peptides, peptide sequences are described in Table 1. The antibody only reacted with the standard peptide and produced a dose-dependent response (Fig. 1a), increased elongated and truncated peptide concentrations only resulted in minor optical density (OD) displacements (Fig. 1b). These data suggest that GPDPLQ₁₂₃₇ is highly specific for the neo-epitope and does not cross-react with the CTX-II neo-epitope, i.e. the Truncated -2 peptide.

Technical validation of the GPDPLO₁₂₃₇ **ELISA**. Technical validation tests were performed to evaluate the performance of the GPDPLQ₁₂₃₇ ELISA. The validation steps and the performance are shown in Table 2. Samples measured linearly in the range 0.8-8.0 ng/ml in dilution recovery tests, samples above 7-8 ng/ml displayed reduced linearity. Samples above 8 ng/ml should be diluted further in assay buffer and measured again. Furthermore, despite obtaining highly reproducible signals in pig, horse, sheep and donkey sera, these signals could not be technically validated by dilution recovery. Hence, sera from these species were only used for determining inter- and intra-assay variation. GPDPLQ₁₂₃₇ was, however, abundant in foetal bovine serum (FBS), where it diluted linearly. Human sera, plasma, and urine from healthy donors or RA and OA patients, as well as rat plasma, did not produce sufficient signals for technical validation.

GPDPLO₁₂₃₇ release from bovine articular cartilage explants. The generation of the GPDPLQ₁₂₃₇ neo-epitope from articular cartilage biopsies undergoing inflammatory degradation was tested in samples derived from bovine cartilage explant (BEX) cultures stimulated with oncostatin M (OSM) and tumour necrosis factor α (TNF α , OT). Although there was substantial variation in GPDPLQ₁₂₃₇ induction by OT in terms

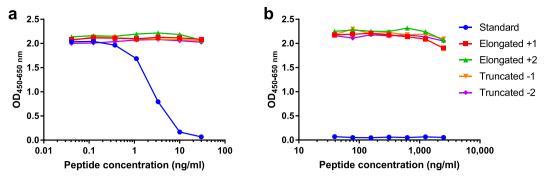


Figure 1. Specificity of the GPDPLQ₁₂₃₇ monoclonal antibody. Monoclonal antibody reactivity towards the standard peptide as well as elongated and truncated standard peptides was tested in the competitive GPDPLQ₁₂₃₇ ELISA. The OD displacement is shown in the peptide concentration range of the standard curve (**a**) as well as at increased peptide concentrations (**b**). Results are shown as $OD_{450-650 \text{ nm}}$ as a function of peptide concentration.

Technical validation procedure	GPDPLQ ₁₂₃₇ performance
Detection range (LLOD-ULOD)	0.7–19.4 ng/ml
Estimated linear range	0.8–8.0 ng/ml
Dilution recovery in linear range †	106% (99–119%)
Inter-assay variation [†]	6% (5-8%)
Intra-assay variation [†]	2% (0-2%)
Freeze/thaw recovery (4 cycles) [†]	100% (99–102%)

Table 2. Technical validation of the GPDPLQ₁₂₃₇ competitive ELISA. [†]Percentages are reported as means with ranges shown in brackets. Abbreviations: LLOD, lower limit of detection; ULOD, upper limit of detection.

of magnitude and kinetics between the three performed experiments, the overall inhibition patterns were the same as in the following representative experiment (Fig. 2). GPDPLQ₁₂₃₇ data from an additional experiment is shown in Supplementary Fig. S1. The GPDPLQ₁₂₃₇ neo-epitope was detected in the OT culture supernatants after approximately 10 days of culture (Fig. 2a), but GPDPLQ₁₂₃₇ levels remained at baseline levels when adding the broad-spectrum MMP inhibitor GM6001 to the culture. The OT response was not inhibited by the cysteine protease inhibitor E-64, albeit with considerable variation between explants, but there was a delay in GPDPLQ₁₂₃₇ release in the E-64 condition relative to the OT condition and the E-64 condition had higher GPDPLQ₁₂₃₇ levels than the OT condition at the final time points, findings that were also present in additional experiments (Supplementary Fig. S1). The apparent decrease in GPDPLQ₁₂₃₇ levels in the OT condition at the final time points compared to earlier time points (Fig. 2a) was, however, not a consistent finding (Supplementary Fig. S1). The GPDPLQ₁₂₃₇ area under the curve (AUC) analysis (Fig. 2c) demonstrates that OT induced a statistically significant increase in GPDPLQ₁₂₃₇ levels (129.9 \pm 58.6 ng/ml, P = 0.034) that was fully inhibited by GM6001 $(10.6 \pm 1.6 \text{ ng/ml}, P = 0.015)$, but not by E-64 (367.2 $\pm 212.2 \text{ ng/ml}, P > 0.999)$. Similar results were also found in additional experiments (Supplementary Fig. S1). CTX-II release, measured at a selection of relevant time points, followed a very similar pattern to that of GPDPLQ₁₂₃₇ release (Fig. 2b); OT induced CTX-II release while GM6001, but not E-64, inhibited the CTX-II induction to baseline levels. Unlike GPDPLQ₁₂₃₇ release, CTX-II release in the E-64 condition did not appear to have a delayed onset relative to the OT condition. Increased CTX-II levels were observed in the E-64 condition at late time points, similarly to the GPDPLQ₁₂₃₇ data. The CTX-II AUC analysis (Fig. 2d) reveals that OT induced a marked increase in CTX-II levels, albeit not statistically significant (5695.0 \pm 1608.0 pg/ml, P = 0.345), that was fully inhibited by GM6001 (283.2 \pm 0.0 pg/ml, P = 0.010), but not by E-64 (23309.0 ± 14381.0 pg/ml, P > 0.999).

GPDPLQ₁₂₃₇ **release from articular cartilage resorption.** The release of GPDPLQ₁₂₃₇ from human osteoclasts cultured on bovine non-calcified articular cartilage was measured in three individual resorption assays. As demonstrated in a representative experiment, the absolute GPDPLQ₁₂₃₇ levels (Fig. 3a) markedly increased in the supernatants of cultures with osteoclasts compared to background levels in cultures without osteoclasts, as did CTX-II levels (Fig. 3b). The background measurements were subtracted from the osteoclast measurements from each assay, and the background-subtracted data were subsequently normalised to the vehicle condition mean. The normalised data from each condition were averaged and pooled with the data from the other resorption assays to generate the final data set, where n corresponds to the number of resorption assays per condition. The background-subtracted GPDPLQ₁₂₃₇ release (Fig. 3c) was inhibited by both E-64 (72.1 ± 14.0% inhibition, P = 0.007) and GM6001 (75.5 ± 13.7% inhibition, P = 0.006) compared to the vehicle condition. The GPDPLQ₁₂₃₇ levels were reduced further by E-64/GM6001 combined (91.5 ± 0.7% inhibition, P = 0.004), although this was not statistically significant compared to E-64 or GM6001 alone (P = 0.659 and P = 0.769).

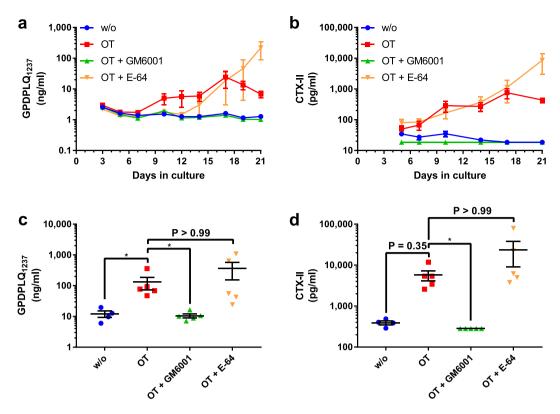


Figure 2. GPDPLQ₁₂₃₇ and CTX-II release into BEX supernatants upon stimulation and inhibition of MMP activity. The BEX cultures were untreated (w/o), treated with OSM and TNF α (OT) without protease inhibitors, or OT with the protease inhibitors GM6001 (OT + GM6001) or E-64 (OT + E-64). GPDPLQ₁₂₃₇ (**a**,**c**) and CTX-II (**b**,**d**) data from one representative experiment are presented as the mean biomarker levels per group at different time points throughout the experiment (**a**,**b**) and as the AUCs of each individual explant and their means (**c**,**d**). Error bars represent the SEMs. Statistical significance relative to the OT condition in the AUC data is indicated by *P < 0.05.

As expected, GPDPLQ₁₂₃₇ was not detectable above background levels in supernatants from resorbed cortical bone (data not shown). CTX-II release (Fig. 3d) in the cartilage resorption cultures was inhibited by GM6001 (87.0 \pm 6.9% inhibition, P = 0.019) but not by E-64. E-64/GM6001 combined did not result in further inhibition of CTX-II release (61.1 \pm 11.6% inhibition, P = 0.246). These findings indicate that the release of GPDPLQ₁₂₃₇ and CTX-II neo-epitopes from articular cartilage undergoing osteoclast-mediated degradation are not mediated by the same proteases; while we demonstrated that both MMPs and cysteine proteases contribute to GPDPLQ₁₂₃₇ release during osteoclastic cartilage resorption, only MMPs contributed to CTX-II release.

Immunolocalization of the GPDPLQ₁₂₃₇ **neo-epitope.** Tibia obtained from young healthy Lewis rats were used to immunolocalise the GPDPLQ₁₂₃₇ neo-epitope. GPDPLQ₁₂₃₇ staining was observed at the articular cartilage and the growth plate, while no staining was observed when using a negative control antibody (Fig. 4). The articular cartilage stained for GPDPLQ₁₂₃₇ in the ECM and in pericellular spaces throughout all cartilage layers, albeit at varying intensity, while the growth plates stained mainly between the columns of proliferating cells—not in the immediate pericellular ECM—and at the proximal bone–cartilage interface. Intermittent growth plate staining was also observed in the ossification zone.

Discussion

There is an unmet need for novel biomarkers of cartilage degradation for rheumatic diseases, and for blood-based biomarkers in particular. The EKGPDPLQ neo-epitope has been previously described in patent literature^{26–28}, but to our knowledge no peer-reviewed articles measuring EKGPDPLQ release from cartilage, calcified or non-calcified, have been published. In patent literature, the EKGPDPLQ neo-epitope has been described to represent the degradation of calcified cartilage by cathepsin K^{26-28} , if analysed in blood. However, using the competitive GPDPLQ₁₂₃₇ ELISA we could quantify both cysteine protease- and MMP-mediated release of the EKGPDPLQ neo-epitope from different models of non-calcified cartilage degradation *in vitro*.

In our study, we detected MMP-derived GPDPLQ₁₂₃₇ release from non-calcified cartilage in a model of inflammatory cartilage degradation. In the BEX OT model, the GPDPLQ₁₂₃₇ protease profile—with inhibition by GM6001 and no inhibition by E-64—and kinetics were similar to those of CTX-II. These CTX-II results also closely resembled previously reported CTX-II^{37,38} and C2M³⁹ data in the BEX OT model. The C2M neo-epitope, like CTX-II, is generated from type II collagen degradation predominantly by MMPs³⁹. C2M has previously been shown to be able to distinguish between healthy subjects and patients or between subsets of patients with

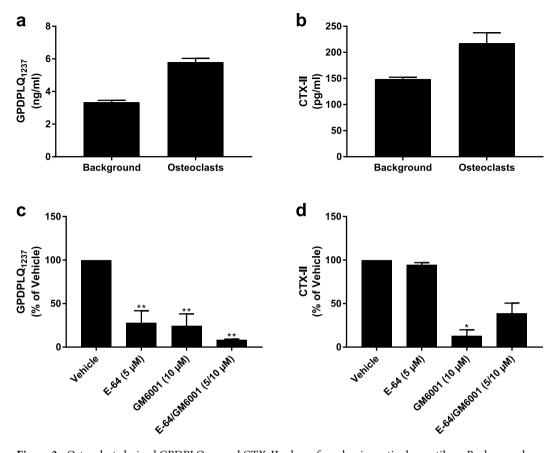


Figure 3. Osteoclast-derived GPDPLQ₁₂₃₇ and CTX-II release from bovine articular cartilage. Background biomarker levels of GPDPLQ₁₂₃₇ (**a**) and CTX-II (**b**) measured in media from cultures with only cartilage (Background, n = 2) and cultures with osteoclasts and cartilage (Osteoclasts, n = 6) in one representative resorption assay, presented as means with SEMs. Background levels were subtracted from the GPDPLQ₁₂₃₇ and CTX-II levels in each individual resorption assay followed by normalisation to the vehicle condition. The normalised data were averaged followed by pooling of the data, resulting in n = 2 for E-64/GM6001 and n = 3 for other conditions where n is the number of resorption assays featuring the condition, and the background-reduced normalised GPDPLQ₁₂₃₇ (**c**) and CTX-II (**d**) levels are presented as the means of all resorption assay means where error bars represent the SEMs. Statistical significance relative to the vehicle condition in the background-reduced data is indicated by *P < 0.05 and **P < 0.01.

arthritic diseases, such as $OA^{39,40}$, RA^{41-43} and ankylosing spondylitis^{44,45}, and is therefore a relevant biomarker for assessing type II collagen resorption in *in vitro* models. Interestingly, CTX-II release is typically increased by E-64 treatment in the BEX OT model^{37,38}, whereas C2M release is not increased to a statistically significant degree³⁹, but we did not detect any significant differences in GPDPLQ₁₂₃₇ or CTX-II release between the OT and the OT+ E-64 conditions. The overall effect of E-64 on both GPDPLQ₁₂₃₇ and CTX-II AUC release was highly variable and not statistically significant, indicating that cysteine proteases did not have a net effect that played any major role in release of either biomarker in this model, although it is possible that specific cysteine proteases may contribute differently. The lack of overall inhibition of GPDPLQ₁₂₃₇ release by E-64 is likely due to low cysteine protease activity in this model, even though cathepsin K induction has previously been demonstrated in the BEX OT model by immunohistochemistry³⁷. Similarly to our findings, the cathepsin K-specific type II collagen neo-epitope C2K77 has been reported to not be detectable above background levels in OT-stimulated equine cartilage explants⁴⁶, indicating that cathepsin K-mediated type II collagen degradation—and potentially the contribution of other cysteine proteases as well—is low in this model.

The reason for the apparent delayed release of GPDPLQ₁₂₃₇ in the E-64 condition as well as the increased GPDPLQ₁₂₃₇ and CTX-II levels at late time points, compared to the OT condition, remains unclear. Parts of these effects are likely explained by the large inter-explant variation that is common for this model system, although we did observe both the delayed GPDPLQ₁₂₃₇ onset and the elevated biomarker levels at final time points in several experiments. The delay in GPDPLQ₁₂₃₇ release did not appear to affect overall biomarker levels over the culture period used in this experiment, and the biological relevance of this putative delay is therefore uncertain. Interestingly, no corresponding delay was observed in CTX-II release upon E-64 treatment, suggesting that there might be a difference between GPDPLQ₁₂₃₇ and CTX-II in the kinetics of biomarker release at intermediate time points that may be related to cysteine protease activity. Speculatively, the delay may be a result of E-64 temporarily inhibiting early cysteine protease-mediated biomarker release, which is detectable by the GPDPLQ₁₂₃₇ assay but

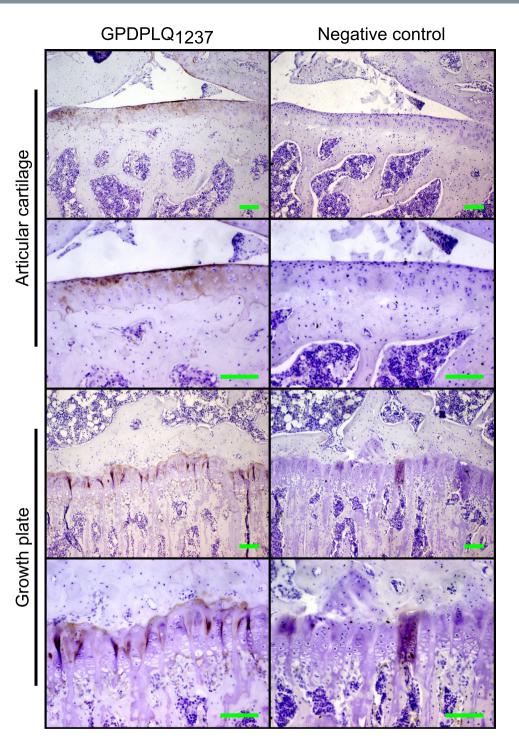


Figure 4. Immunolocalization of the GPDPLQ₁₂₃₇ neo-epitope in rat tibia. Representative micrographs of rat tibial articular cartilage and growth plate stained with the NB427-5G11-4W3 GPDPLQ₁₂₃₇ (GPDPLQ₁₂₃₇) monoclonal antibody (brown) or a negative control antibody, counterstained with Mayer's haematoxylin (purple). Micrographs were obtained using a 10x (first and third row) and a 20x (second and fourth row) objective. Scale bars represent 100 μ m.

not by the CTX-II assay, but the inhibition ultimately fails and results in rapid onset of cysteine protease-mediated cartilage degradation and GPDPLQ_{1237} release.

The increased GPDPLQ₁₂₃₇ release at late time points in the E-64 condition may potentially be mediated by inhibition of cathepsin B. Cathepsin B is a dipeptidyl carboxypeptidase⁴⁷ that is upregulated in interleukin-1 (IL-1)⁴⁸ and TNF α /IL-1 β -stimulated⁴⁹ chondrocytes and may potentially process the EKGPDPLQ neo-epitope into the EKGPDP neo-epitope. Although cathepsin B is known to generate CTX-II during *in vitro* cleavage—while cathepsins D, K, L and S do not—and may even be able to degrade CTX-II neo-epitopes to some extent⁵⁰, it is unclear how much cathepsin B actually contributes to CTX-II release in the BEX OT model. As demonstrated in

our study, cysteine protease inhibition with E-64 does not reduce CTX-II release in the BEX OT model but rather appears to have no effect or to increase it, as has been previously reported^{37,38}, indicating that the contribution of cathepsin B to type II collagen degradation, CTX-II release, GPDPLQ₁₂₃₇ degradation, and/or CTX-II degradation in this model is modest compared to the net neutral or positive effects that cysteine protease inhibition appears to have on CTX-II levels at late time points. Furthermore, the overall lack of effect of E-64 in GPDPLQ₁₂₃₇ and CTX-II AUC data and the increased release of both GPDPLQ1237 and CTX-II at late time points suggest that the late increase of GPDPLQ₁₂₃₇ levels is not a result of inhibited conversion of the GPDPLQ₁₂₃₇ epitope into CTX-II by cysteine proteases—as E-64 treatment only resulted in delayed GPDPLQ1237 induction while it did not increase and decrease overall GPDPLQ₁₂₃₇ and CTX-II release, respectively. While cysteine proteases do not appear to have a significant net effect on overall GPDPLQ₁₂₃₇ and CTX-II release, cathepsin B may still convert GPDPLQ₁₂₃₇ into CTX-II in this model. If cathepsin B mediates both GPDPLQ₁₂₃₇ and CTX-II degradation at similar rates it is possible that neither of these effects would be detectable during cysteine protease inhibition, or indeed cathepsin B inhibition. The potential role of cathepsin B in conversion of GPDPLQ₁₂₃₇ into CTX-II is interesting and could be investigated further using cathepsin B-specific inhibitors such as CA074⁵¹, potentially in BEX cultures stimulated with IL-1, but this was considered to be outside the scope of the current study as cysteine protease inhibition did not appear to have substantial impact on the overall release of either biomarker in the BEX OT model.

As for the observed decline in biomarker release in the OT condition at very late time points, this is similar to previously reported results in other type II collagen neo-epitopes at very late time points, both in OT and OT+E-64 conditions^{38,39}, but it was not a consistent finding in our study, underscoring the somewhat variable nature of cartilage explant cultures. As it was not consistent and was observed with both GPDPLQ₁₂₃₇ and CTX-II, albeit to a lesser extent with CTX-II, the decline is unlikely to be a result of interconversion of the neo-epitopes by cathepsin B or to represent important changes to type II collagen degradation over time. The decline may rather be related to e.g. destruction of the cartilage or depletion of the ECM, as has been previously described in this model³⁷, or the substantial inter-explant variation that is often observed in the BEX OT model.

In our model of osteoclastic cartilage resorption, GPDPLQ₁₂₃₇ release was found to be mediated by both cysteine proteases and MMPs. Cathepsin K is the predominant cathepsin in osteoclasts^{52,53} and it has also been demonstrated to be crucial for type I collagen degradation by osteoclasts during normal bone resorption⁵²⁻⁵⁶. Furthermore, the cathepsin K-specific antiresorptive odanacatib inhibits the majority of CTX-I release (approximately 80% inhibition) during bone resorption in vitro^{55,56} with an effect size similar to that of E-64⁵⁷. Osteoclasts are also thought to utilize several other cathepsins and cysteine proteases for bone resorption^{58,59}, but cathepsin K is clearly essential for normal resorptive function and type I collagen resorption even though MMPs can compensate to some extent when cathepsin K is inhibited⁶⁰. Accordingly, it is likely that the cysteine protease-mediated GPDPLQ₁₂₃₇ release from cartilage is mediated at least in part by cathepsin K, but the precise cysteine proteases underlying GPDPLQ₁₂₃₇ release during cartilage resorption would need to be verified using specific inhibitors, like odanacatib, or cathepsin K-deficient osteoclasts. Contrary to GPDPLQ₁₂₃₇ release, CTX-II release from cartilage resorption was mediated by MMPs but not by cysteine proteases. This is similar to previous studies of CTX-II release from cartilage digested by MMPs or cathepsins D, K, L and S in vitro^{50,61}, although CTX-II has been shown to be generated by cathepsin B⁵⁰. Thus, GPDPLQ₁₂₃₇ appears to be a biomarker that is generated through direct osteoclastic cartilage resorption mediated by cysteine proteases—likely in part by cathepsin K and MMPs whereas CTX-II release is primarily mediated by MMPs in this model, similarly to previous studies of C2M release in this model³⁶.

The potential contribution of MMPs, cysteine proteases and articular cartilage to GPDPLQ₁₂₃₇ neo-epitope levels may differ in vivo. GPDPLQ₁₂₃₇ levels in human serum, plasma and urine, obtained from healthy donors, OA patients, or RA patients, or rat plasma were too low to perform technical validation of the assay in these matrices. Further work would be required to validate the neo-epitope's relevance in human samples. Although we do not know the reason underlying the lack of substantial GPDPLQ₁₂₃₇ signals in serum, plasma, and urine in this study, it may be due to insufficient sensitivity of the ELISA or processing of the EKGPDPLQ neo-epitope into the EKGPDP neo-epitope²⁶⁻²⁸, followed by excretion, resulting in low levels or transient presence of the EKGPDPLQ neo-epitope in serum and plasma. This would require further testing and validation, potentially on more sensitive assay platforms or using a redesigned assay that may allow for improved sensitivity towards the analyte. However, preliminary attempts to change the assay buffer or assay design, using e.g. a sandwich ELISA instead of a competitive ELISA, did not improve sensitivity towards the analyte in serum or plasma. Alternatively, GPDPLQ₁₂₃₇ may be measured in SF which may be a more biologically relevant matrix, although the patents proposed otherwise²⁶⁻²⁸. SF measurements were outside the scope of the current study as we aimed to develop a blood- or urine-based biomarker, due to the limited availability of SF and its invasive sampling method, and GPDPLQ₁₂₃₇ presence in SF would need to be investigated further. Other elongated or truncated CTX-II neo-epitopes²⁶⁻²⁸ may also be relevant targets in the search for CTX-II-like neo-epitopes in blood.

Despite the lack of GPDPLQ₁₂₃₇ in human and rat body fluids, our immunohistochemical analysis of GPDPLQ₁₂₃₇ presence in rat tibia suggests that the neo-epitope can be generated through naturally occurring cartilage turnover, at least in the joints of young adult rats. GPDPLQ₁₂₃₇ was observed in the articular cartilage ECM in healthy Lewis rats. CTX-II is primarily present in articular cartilage during disease, e.g. collagen-induced arthritis, in rats^{18,62} but it has also been described in some healthy rats⁶². In our study, GPDPLQ₁₂₃₇ was present at various levels in all layers of the articular cartilage ECM in healthy rats, unlike CTX-II which is only present in some healthy rats⁶², which suggests that articular cartilage GPDPLQ₁₂₃₇ may be a feature of physiological cartilage turnover in rats, whereas CTX-II is a better indicator of pathology. Additionally, the articular cartilage GPDPLQ₁₂₃₇ staining patterns differ from those reported for CTX-II; CTX-II stains the superficial matrix in the upper zone as well as around round and flat chondrocytes in the upper and deep zones⁶², although it is unclear which treatments the rats that were used in the reported illustrations were subjected to in this report.

Similar CTX-II staining patterns have also been reported in human OA knees²⁹. Interestingly, CTX-II frequently stains the bone–cartilage interface in these tissues²⁹, but human growth plate stainings were not available for comparison to our GPDPLQ₁₂₃₇ findings. When examining the rat growth plates, we found that GPDPLQ₁₂₃₇ primarily stained the proximal bone–cartilage interface and the ECM between the columns of proliferating chondrocytes, but not the pericellular matrix. We also observed intermittent GPDPLQ₁₂₃₇ staining in the ossification zone, where osteoclasts degrade the calcified cartilage during endochondral ossification⁶³. Growth plates in both healthy rats and arthritic rats can stain positive for the CTX-II neo-epitope, but they do so in a different pattern than GPDPLQ₁₂₃₇⁶². CTX-II staining was reported primarily in the pericellular ECM around the chondrocytes in the illustrative growth plate staining⁶². Although the GPDPLQ₁₂₃₇ and CTX-II staining patterns clearly differ, the presence of both GPDPLQ₁₂₃₇ and CTX-II at healthy rat growth plate suggests that they may both be generated as part of endochondral ossification during bone growth. This conclusion is further supported by both GPDPLQ₁₂₃₇ and CTX-II¹⁷ being abundant in FBS, as endochondral ossification is an important part of bone development during foetal development. Similarly, both CTX-II and C2M are high in serum from one-month-old rats but decrease significantly over time as they age⁶⁴.

Other type II collagen neo-epitope antibodies, such as C2M and C2K77 have also been described to be present in OA cartilage. C2M is known to be present in both non-calcified and calcified cartilage beneath cartilage surface irregularities, surface erosions, and cartilage lesions³⁹. Similarly, C2K77 staining is increased in equine OA cartilage sections with focal cartilage degeneration and partial thickness fissures⁴⁶, as well as in areas of proteoglycan loss⁶⁵. C2K77 staining has also been demonstrated in human knee OA cartilage lesions where it mostly stained the matrix around chondrocytes in the lesional area of the tissue, whereas little staining was observed in non-lesion areas or with a negative control antibody⁶⁵. These immunohistochemical stainings of the type II collagen neo-epitopes CTX-II, C2M and C2K77 demonstrate how they can be used to further our understanding of disease processes in different parts of joints, but descriptions of C2M, C2K77 and CTX-II stainings of healthy cartilage are limited. In our study we observed substantial GPDPLQ₁₂₃₇ staining in both cartilage and growth plates in healthy rat specimen. However, our immunohistochemical analysis does demonstrate that GPDPLQ₁₂₃₇ is generated biologically and is present in rat tissues, even though we were not able to measure and validate GPDPLQ₁₂₃₇ is rat plasma. Immunohistochemical studies using human tissues would be of interest to investigate if GPDPLQ₁₂₃₇ staining is also a feature of mature human cartilage or if it is specific to diseased human cartilage.

In conclusion, the results of our study indicate that the GPDPLQ₁₂₃₇ neo-epitope—EKGPDPLQ—is a non-calcified cartilage degradation biomarker that can be generated during osteoclastic resorption and inflammatory degradation by cysteine proteases and/or MMPs *in vitro*, and during physiological cartilage degradation in rats. GPDPLQ₁₂₃₇ is a promising biomarker for studies of cartilage degradation *in vitro*, but additional work is required to validate the relevance of the neo-epitope *in vivo*.

Methods

Target selection and peptide sequences. The GPDPLQ₁₂₃₇ target sequence, as described in patents^{26–28}, is a C-terminal neo-epitope formed by proteolytic cleavage of type II collagen: $_{1230}$ EKGPDPLQ \downarrow YMR₁₂₄₀. A sequence of 8 amino acids adjacent to the cleavage site ($_{1230}$ EKGPDPLQ \downarrow_{1237}) was chosen as the target. The sequence was analysed for homology to other human and animal proteins using NPS@: Network Protein Sequence Analysis PattInProt search⁶⁶ with the UniProtKB/Swiss-prot database.

Synthetic peptides (GenScript, Piscataway, NJ) were used for monoclonal antibody production and validation of the ELISA. An overview of the peptides used is shown in Table 1. A biotinylated peptide (biotin-EKGPDPLQ) was used as a coating peptide for streptavidin-coated microtiter plates. Standard curves were generated using the standard peptide (EKGPDPLQ). The specificity of the antibody was tested using two elongated and two truncated standard peptides, that have been previously described^{26–28}. An immunogenic peptide was generated by covalent cross-linkage of a standard peptide with an added N-terminal cysteine residue to keyhole limpet haemocyanin (KLH) carrier protein using succinimidyl 4-(*N*-maleimidomethyl)cyclohexane-1-carboxylate (SMCC, Thermo Scientific, Waltham, MA).

Generation of monoclones. Four to six-week-old female Balb/C mice were housed (five mice per cage, standard wood chips bedding enriched with huts, nesting material and sticks) at the Nordic Bioscience (Beijing) animal facility (21-23 °C, 55–65% relative humidity, 12 hours light/dark cycle) with *ad libitum* access to food and water. The mice were immunised by subcutaneous injection of 200 µl emulsified antigen containing 60 µg immunogenic peptide mixed with Freund's complete adjuvant (Sigma-Aldrich, St. Louis, MO). Repeated immunizations with 30 µg immunogenic peptide mixed with Freund's incomplete adjuvant (Sigma-Aldrich, St. Louis, MO) were performed in 2-week intervals until stable serum titres were obtained. The mice with the highest titre and highest reactivity towards the target, tested in competitive ELISAs using the biotinylated, standard, elongated and truncated peptides, were allowed to recover for four weeks followed by a booster injection of 50 µg immunogenic peptide in 100 µl 0.9% NaCl solution, administered intraperitoneally. When serum titres were high and stable, the mice were euthanised by concussion followed by cervical dislocation and isolation of their spleens. Hybridoma cells were generated through fusion of isolated splenocytes with SP2/0 myeloma cells as previously described⁶⁷. Limited dilution procedures were used to generate monoclonal cultures.

Specificity of monoclonal antibodies. The reactivity of the monoclonal supernatants was evaluated in a competitive ELISA by testing for OD displacement using the standard peptide and lack of displacement by elongated and truncated standard peptides, using 4 ng/ml biotinylated peptide on streptavidin-coated microtiter plates (Roche, Basel, Switzerland, cat. 11940279). The clones with the best OD displacement profiles (data not

shown) were used to generate antibodies that were purified using protein G columns according to the manufacturer's instructions (GE Healthcare Life Sciences, Little Chalfont, UK, cat. 17-0404-01).

GPDPLO₁₂₃₇ **assay protocol.** Optimal assay buffer, incubation times, temperature, and concentrations of antibody and peptides were determined and the finalised GPDPLQ₁₂₃₇ competitive ELISA protocol is described below. A 96-well streptavidin plate was coated with 100 µl of 625 pg/ml biotinylated peptide dissolved in coating buffer (10 mM phosphate-buffered saline [PBS] with 1% bovine serum albumin and 0.1% Tween 20 [PBS-BTB], 8 g/l NaCl, pH 7.4) and incubated for 30 minutes at 20 °C. 20 µl of peptide calibrators prepared as a 3-fold dilution series of 30.03 ng/ml standard peptide, controls or samples were added to the appropriate wells followed by $100\,\mu$ l of antibody NB427-5G11-4W3 dissolved in assay buffer (50 mM PBS-BTB, 8 g/l NaCl, pH 7.4) to a concentration of 526 ng/ml and incubation for 20±1 hours at 4 °C. 100 µl of horseradish peroxidase (HRP)-conjugated rabbit anti-mouse secondary antibody (Jackson ImmunoResearch Laboratories, West Grove, PA, cat. 315-035-045) dissolved in assay buffer to a concentration of 160 ng/ml was added to the wells and incubated for 1 hour at 20 °C. All the above incubation steps included shaking at 300 rpm and were followed by five washes in washing buffer (25 mM Tris, 50 mM NaCl, 0.1% Tween 20, pH 7.2). Finally, 100 µl 3,3',5,5'-tetramethylbenzidine (TMB, Kem-En-Tec Diagnostics, Taastrup, Denmark, cat. 4380) was added and incubated for 15 minutes at 20 °C in the dark with shaking at 300 rpm followed by the addition of $100\,\mu$ l of $0.18\,M\,H_2SO_4$ to stop the reaction. The absorbance was measured at 450 nm with 650 nm as a reference on a VersaMax (Molecular Devices, Sunnyvale, CA) and a standard curve was plotted using a 4-parametric mathematical fit model in SoftMax Pro v.6.3 for Windows (Molecular Devices, Sunnyvale, CA). Cartilage resorption supernatants were diluted 1:1-1:5 as needed, BEX supernatants were diluted approximately 1:1-1:4000 as needed, serum, plasma, and urine samples were measured undiluted. Any samples measuring below the lower limit of detection (LLOD) were assigned the LLOD as their result.

Technical validation. The technical performance of the ELISA was validated in the following tests: LLOD, upper limit of detection (ULOD), measuring range (LLOD-ULOD), linear range and dilution recovery, interand intra-assay variation, specificity, and analyte freeze-thaw stability. The LLOD was determined from three independent runs of 68 samples of the zero standard (i.e. assay buffer). The LLOD was calculated within each run as the mean + three standard deviations and the final LLOD was calculated as the mean of the three runs. The ULOD was determined from five independent runs of the highest standard, measured in duplicates. The mean concentrations of the duplicates, as calculated from the standard curve, were used to calculate a ULOD defined as the mean of the five measurements - three standard deviations. The linear range of the assay was estimated from dilution recovery tests by performing 2-fold dilution series of a heat-inactivated FBS sample and a culture supernatant from a BEX OT culture and calculating the percentage dilution recovery against different starting dilutions. The linear range was considered to be the longest range of measured concentrations that, when adjusted for dilution, had a dilution recovery within $100 \pm 20\%$ of the adjusted concentration of the starting dilution. Interand intra-assay variation was determined from five independents runs of six positive control samples-diluted FBS, diluted BEX OT, and sera from pig, horse, sheep and donkey measured in duplicates—and the inter- and intra-assay coefficients of variation (CVs) of each sample were calculated. Specificity was tested using elongated and truncated standard peptides (Table 1). The analyte stability was determined for one FBS sample and one BEX OT sample for up to four freeze-thaw cycles.

BEX cultures. BEXs were harvested by dissecting the outermost layer of articular cartilage from bovine knee joints and were cultured as previously described⁶⁸. In brief, cartilage explants were placed in 96-well plates and incubated at 37 °C with 5% CO₂ under serum-free conditions (n = 4–5 per condition). Each explant was cultured in 200 µl DMEM/F-12 medium for 21 days, with medium changes every 3–4 days, under one of the following conditions: without cytokines, with 10 ng/ml OSM and 20 ng/ml TNF α (OT) to stimulate MMP activity, OT supplemented with 10 µM of the broad-spectrum MMP inhibitor GM6001⁶⁹, or OT supplemented with 10 µM of the cysteine protease inhibitor E-64⁷⁰. The selection of inhibitors and dosages was made based on their previous use in cartilage explant^{37–39} and cartilage resorption cultures³⁶ as well as their well-established effects on osteoclastic bone resorption⁵⁷.

Articular cartilage resorption assays. Osteoclast-mediated degradation of bovine articular cartilage was performed as previously described³⁶. In brief, bovine articular cartilage was isolated from the femoral condyles of three bovine knees, using a biopsy punch and a scalpel, followed by immersion in liquid nitrogen to render the cartilage metabolically inactive and subsequent fixation of the tissue in 70% ethanol at 4°C for at least 5 days prior to use. Mature human osteoclasts were generated from CD14⁺ monocytes isolated from the peripheral blood of anonymised blood donors as previously described^{71,72}. Osteoclasts were seeded on the cartilage in 96-well culture plates at a density of 1.0×10^5 cells/well and were cultured at 37 °C and 5% CO₂ in alpha minimum essential medium (α MEM) containing 10% heat-inactivated FBS, 100 units/ml penicillin, 100 µg/ml streptomycin, 388 µg/l thymidine, 25 ng/ml macrophage colony-stimulating factor (M-CSF, R&D Systems, Minneapolis, MN) and 25 ng/ ml receptor activator of nuclear factor kappa-B ligand (RANKL, R&D Systems). The medium was supplemented with 5µM E-64, 10µM GM6001, 5µM E-64 and 10µM GM6001 combined (E-64/GM6001), or a DMSO vehicle (1:2000 in medium), all from Sigma-Aldrich (St. Louis, MO). The selection of inhibitors and dosages was made based on their well-established effects on osteoclastic bone resorption⁵⁷ as well as their previous use in cartilage resorption³⁶ and cartilage explant³⁷⁻³⁹ cultures. Cartilage without osteoclasts was cultured as background samples for biomarker measurements. The media were changed by demi-depletion after 24 hours. After an additional 3-4 days the media were collected and stored at -20 °C until analysis. Three independent resorption assays (n = 3-6 per condition) were performed in this study.

CTX-II. CTX-II released from resorbed cartilage and BEX cultures was measured in the culture supernatants using the Serum Pre-Clinical CartiLaps ELISA (IDS, The Boldons, UK), according to the manufacturer's instructions. Cartilage resorption supernatants were diluted 1:1–1:8. BEX supernatants were measured for CTX-II at a selection of relevant time points, based on the GPDPLQ₁₂₃₇ induction data (days 5, 7, 10, 14, 17, and 21), with samples diluted 1:5–1:500 as needed. Any samples measuring below the detection limit, as specified by the manufacturer, were assigned the detection limit (3.7 pg/ml) as their result. The detection limit was adjusted to 18.5 pg/ml for the BEX samples to compensate for the lowest dilution factor used being 1:5.

Immunohistochemistry. Knees from young adult (19–21 weeks) female Lewis rats were fixed in 10% formalin for 2 weeks followed by decalcification in 15% EDTA for 7 weeks. The decalcified knees were infiltrated with paraffin using a Tissue-Tek VIP 5 Jr. (Sakura Finetek, Alphen aan den Rijn, The Netherlands), embedded in paraffin, and cut into 5-6 µm thick coronal sections using a HM 360 microtome (Microm International GmbH, Walldorf, Germany). The sections were deparaffinised and subjected to heat-induced epitope retrieval in sodium citrate buffer (10 mM tri-sodium citrate dihydrate, 0.05% Tween 20, pH 6) at 60 °C overnight. The following day the sections were washed in tris-buffered saline (TBS) with 0.1% Tween 20 and 1% Triton X-100 (TBS-TT, used for all wash steps), blocked for endogenous peroxidase activity using 1.2% hydrogen peroxide in 70% ethanol, washed, incubated with blocking solution (0.5% casein sodium salt in TBS-TT) for 20 minutes at room temperature, and incubated with the NB427-5G11-4W3 monoclonal antibody or a negative control mouse antibody (Agilent Technologies, Santa Clara, CA, cat. X0931)-dissolved in blocking solution to 5.1 µg/ml and 15.6 µg/ ml respectively—overnight at 4 °C. Then the sections were washed, incubated with an HRP-conjugated secondary antibody (Agilent Technologies, Santa Clara, CA, cat. K4001) for 30 minutes at room temperature, washed, incubated with diaminobenzidine (DAB) chromogen, rinsed in tap water, counterstained with Mayer's haematoxylin, rinsed and mounted. Digital micrographs of the tibial articular cartilage and growth plate were obtained with an Olympus DP71 digital camera mounted on a BX-60 microscope with 10x and 20x objectives using the Olympus cellSens software (Olympus, Center Valley, PA). The brightness and contrast of each individual image was adjusted for clarity using the Fiji⁷³ Brightness/Contrast tool in the automatic mode.

Statistical analysis. The BEX biomarker data were used to calculate the AUC of each explant. The group AUCs were analysed using the Kruskal-Wallis test with Dunn's multiple comparisons test against the OT condition. Statistical significance was considered to be P < 0.05. Means \pm standard errors of the mean (SEMs) and P-values from the multiple comparisons tests are reported in the results section. Statistical significance is indicated with symbols in the figure where *P < 0.05. Statistical comparisons and plotting of graphs were performed in GraphPad Prism v7.01 for Windows (GraphPad Software, La Jolla, CA). Graphs represent group means and SEMs.

For the statistical analysis of cartilage resorption assays, the background biomarker measurements were subtracted from all other conditions within the same resorption assay, n = 3-6 per condition. The background-subtracted data were normalised to the vehicle condition mean. The mean percentages of each condition from each resorption assay were pooled for analysis (n = 2 for E-64/GM6001 and n = 3 for other conditions, n represents the number of resorption assays the condition was used in), reported as percentages of the vehicle \pm SEMs. The data were analysed with a two-tailed one-way analysis of variance (ANOVA), with the assumption that the data were normally distributed, followed by either Tukey (for GPDPLQ₁₂₃₇) or Dunnett's T3 (for CTX-II) post hoc tests depending on whether or not equal variances could be assumed based on Levene's test. Statistical significance was considered to be P < 0.05. Significance levels from the post hoc tests, relative to the vehicle condition (unless otherwise specified), are reported with P-values in the results section and using symbols in the figure where *P < 0.05 and **P < 0.01. Statistical comparisons were performed using IBM SPSS Statistics v24.0.0.0 64-bit edition for Windows (IBM, Armonk, NY). Graphs were plotted in GraphPad Prism v7.01 for Windows (GraphPad Software, La Jolla, CA) and represent group means and SEMs.

Ethical approval. Mouse immunisation procedures were approved by the Beijing Administration Office of Laboratory Animal and the animal ethics committee of Nordic Bioscience (Beijing).

The use of peripheral blood from anonymised blood donors, obtained from a blood bank, for osteoclast cultures in this study was covered by the general ethical approval for research use of donor material in accordance with the Transfusion Medicine Standards (TMS) of the Danish Society of Clinical Immunology (DSKI). Informed consent to participate as an anonymised healthy control was obtained from each of the participating donors in writing as part of standard practices at the blood bank at Rigshospitalet, Copenhagen University Hospital, Copenhagen, Denmark. Donors participated in compliance with the Helsinki Declaration.

Human serum, plasma, and urine from OA patients, RA patients, and healthy donors for GPDPLQ₁₂₃₇ measurements were commercially sourced.

Data Availability

The datasets generated during and/or analysed during the current study are available from the corresponding author on reasonable request.

References

- 1. Valdes, A. M. *et al.* Large scale meta-analysis of urinary C-terminal telopeptide, serum cartilage oligomeric protein and matrix metalloprotease degraded type II collagen and their role in prevalence, incidence and progression of osteoarthritis. *Osteoarthr. Cartil.* **22**, 683–9 (2014).
- Lohmander, L. S., Atley, L. M., Pietka, T. A. & Eyre, D. R. The release of crosslinked peptides from type II collagen into human synovial fluid is increased soon after joint injury and in osteoarthritis. *Arthritis Rheum.* 48, 3130–9 (2003).
- Leung, Y. Y., Huebner, J. L., Haaland, B., Wong, S. B. S. & Kraus, V. B. Synovial fluid pro-inflammatory profile differs according to the characteristics of knee pain. Osteoarthr. Cartil. 25, 1420–1427 (2017).

- 4. Krabben, A., Huizinga, T. W. & Mil, A. H. Biomarkers for radiographic progression in rheumatoid arthritis. *Curr. Pharm. Des.* 21, 147–69 (2015).
- 5. Reijman, M. et al. A new marker for osteoarthritis: cross-sectional and longitudinal approach. Arthritis Rheum. 50, 2471-8 (2004).
- Kraus, V. B. et al. Predictive validity of biochemical biomarkers in knee osteoarthritis: data from the FNIH OA Biomarkers Consortium. Ann. Rheum. Dis. 76, 186–195 (2017).
- 7. Garnero, P. et al. Cross sectional evaluation of biochemical markers of bone, cartilage, and synovial tissue metabolism in patients with knee osteoarthritis: relations with disease activity and joint damage. Ann. Rheum. Dis. **60**, 619–26 (2001).
- Jordan, K. M. et al. Urinary CTX-II and glucosyl-galactosyl-pyridinoline are associated with the presence and severity of radiographic knee osteoarthritis in men. Ann. Rheum. Dis. 65, 871–7 (2006).
- Garnero, P., Peterfy, C., Zaim, S. & Schoenharting, M. Bone marrow abnormalities on magnetic resonance imaging are associated with type II collagen degradation in knee osteoarthritis: a three-month longitudinal study. *Arthritis Rheum.* 52, 2822–9 (2005).
- Van Spil, W. E. et al. The ability of systemic biochemical markers to reflect presence, incidence, and progression of early-stage radiographic knee and hip osteoarthritis: data from CHECK. Osteoarthr. Cartil. 23, 1388–97 (2015).
- 11. Thudium, C. S., Löfvall, H., Karsdal, M. A., Bay-Jensen, A.-C. & Bihlet, A. R. Protein biomarkers associated with pain mechanisms in osteoarthritis. J. Proteomics 190, 55–66 (2019).
- Garnero, P. et al. Association of baseline levels of markers of bone and cartilage degradation with long-term progression of joint damage in patients with early rheumatoid arthritis: the COBRA study. Arthritis Rheum. 46, 2847–56 (2002).
- 13. Landewé, R. *et al.* Markers for type II collagen breakdown predict the effect of disease-modifying treatment on long-term radiographic progression in patients with rheumatoid arthritis. *Arthritis Rheum.* **50**, 1390–9 (2004).
- 14. Young-Min, S. et al. Biomarkers predict radiographic progression in early rheumatoid arthritis and perform well compared with traditional markers. Arthritis Rheum. 56, 3236–47 (2007).
- Marotte, H., Gineyts, E., Miossec, P. & Delmas, P. D. Effects of infliximab therapy on biological markers of synovium activity and cartilage breakdown in patients with rheumatoid arthritis. *Ann. Rheum. Dis.* 68, 1197–200 (2009).
- Freeston, J. E. et al. Urinary type II collagen C-terminal peptide is associated with synovitis and predicts structural bone loss in very early inflammatory arthritis. Ann. Rheum. Dis. 70, 331–3 (2011).
- Oestergaard, S. *et al.* The utility of measuring C-terminal telopeptides of collagen type II (CTX-II) in serum and synovial fluid samples for estimation of articular cartilage status in experimental models of destructive joint diseases. *Osteoarthr. Cartil.* 14, 670–9 (2006).
- Oestergaard, S. *et al.* Early elevation in circulating levels of C-telopeptides of type II collagen predicts structural damage in articular cartilage in the rodent model of collagen-induced arthritis. *Arthritis Rheum.* 54, 2886–90 (2006).
- 19. Oestergaard, S. *et al.* Effects of ovariectomy and estrogen therapy on type II collagen degradation and structural integrity of articular cartilage in rats: implications of the time of initiation. *Arthritis Rheum.* **54**, 2441–51 (2006).
- Nielsen, R. H., Christiansen, C., Stolina, M. & Karsdal, M. A. Oestrogen exhibits type II collagen protective effects and attenuates collagen-induced arthritis in rats. *Clin. Exp. Immunol.* 152, 21–7 (2008).
- Duclos, M. E. et al. Significance of the serum CTX-II level in an osteoarthritis animal model: a 5-month longitudinal study. Osteoarthr. Cartil. 18, 1467–76 (2010).
- 22. Kraus, V. B. et al. Urea as a passive transport marker for arthritis biomarker studies. Arthritis Rheum. 46, 420-7 (2002).
- Karsdal, M. A. *et al.* Biochemical markers identify influences on bone and cartilage degradation in osteoarthritis-the effect of sex, Kellgren-Lawrence (KL) score, body mass index (BMI), oral salmon calcitonin (sCT) treatment and diurnal variation. *BMC Musculoskelet. Disord.* 11, 125 (2010).
- 24. Schaller, S. et al. In vitro, ex vivo, and in vivo methodological approaches for studying therapeutic targets of osteoporosis and degenerative joint diseases: how biomarkers can assist? Assay Drug Dev. Technol. 3, 553-80 (2005).
- 25. Christgau, S. et al. Collagen type II C-telopeptide fragments as an index of cartilage degradation. Bone 29, 209-15 (2001).
- 26. Eyre, D. R. Cartilage resorption assays. US6255056B1 (2001).
- 27. Eyre, D. R. Cartilage resorption assays measuring type II collagen fragments. US6348320B1 (2002).
- 28. Eyre, D. R. Synthetic peptides of type II collagen for cartilage resorption assays. US6566492B2 (2003).
- 29. Bay-Jensen, A.-C. *et al.* Biochemical markers of type II collagen breakdown and synthesis are positioned at specific sites in human osteoarthritic knee cartilage. *Osteoarthr. Cartil.* **16**, 615–23 (2008).
- 30. Blaney Davidson, E. N., Vitters, E. L., van der Kraan, P. M. & van den Berg, W. B. Expression of transforming growth factor-β (TGFβ) and the TGFβ signalling molecule SMAD-2P in spontaneous and instability-induced osteoarthritis: role in cartilage degradation, chondrogenesis and osteophyte formation. Ann. Rheum. Dis. 65, 1414–21 (2006).
- 31. Diarra, D. et al. Dickkopf-1 is a master regulator of joint remodeling. Nat. Med. 13, 156-63 (2007).
- Schett, G., Hayer, S., Zwerina, J., Redlich, K. & Smolen, J. S. Mechanisms of Disease: the link between RANKL and arthritic bone disease. Nat. Clin. Pract. Rheumatol. 1, 47–54 (2005).
- 33. Schett, G. et al. Analysis of the kinetics of osteoclastogenesis in arthritic rats. Arthritis Rheum. 52, 3192-201 (2005).
- Stolina, M. et al. RANKL is a marker and mediator of local and systemic bone loss in two rat models of inflammatory arthritis. J. Bone Miner. Res. 20, 1756–65 (2005).
- Goldring, S. R. & Goldring, M. B. Changes in the osteochondral unit during osteoarthritis: structure, function and cartilage-bone crosstalk. Nat. Rev. Rheumatol. 12, 632–644 (2016).
- 36. Löfvall, H. *et al.* Osteoclasts degrade bone and cartilage knee joint compartments through different resorption processes. *Arthritis Res. Ther.* **20**, 67 (2018).
- 37. Sondergaard, B. C. et al. Relative contribution of matrix metalloprotease and cysteine protease activities to cytokine-stimulated articular cartilage degradation. Osteoarthr. Cartil. 14, 738–48 (2006).
- He, Y. et al. The effect of protease inhibitors on the induction of osteoarthritis-related biomarkers in bovine full-depth cartilage explants. PLoS One 10, e0122700 (2015).
- Bay-Jensen, A.-C. et al. Enzyme-linked immunosorbent assay (ELISAs) for metalloproteinase derived type II collagen neoepitope, CIIM—increased serum CIIM in subjects with severe radiographic osteoarthritis. Clin. Biochem. 44, 423–9 (2011).
- 40. Siebuhr, A. S. *et al.* Identification and characterisation of osteoarthritis patients with inflammation derived tissue turnover. *Osteoarthr. Cartil.* **22**, 44–50 (2014).
- Maijer, K. I. et al. Neo-Epitopes—Fragments of Cartilage and Connective Tissue Degradation in Early Rheumatoid Arthritis and Unclassified Arthritis. PLoS One 11, e0149329 (2016).
- 42. Bay-Jensen, A. C. et al. Effect of tocilizumab combined with methotrexate on circulating biomarkers of synovium, cartilage, and bone in the LITHE study. Semin. Arthritis Rheum. 43, 470–8 (2014).
- 43. Bay-Jensen, A. C. *et al.* Early changes in blood-based joint tissue destruction biomarkers are predictive of response to tocilizumab in the LITHE study. *Arthritis Res. Ther.* **18**, 13 (2016).
- Bay-Jensen, A. C. et al. Ankylosing spondylitis is characterized by an increased turnover of several different metalloproteinasederived collagen species: A cross-sectional study. Rheumatol. Int. 32, 3565–3572 (2012).
- 45. Bay-Jensen, A. C. et al. Circulating protein fragments of cartilage and connective tissue degradation are diagnostic and prognostic markers of rheumatoid arthritis and ankylosing spondylitis. PLoS One 8, e54504 (2013).
- Noé, B. et al. C2K77 ELISA detects cleavage of type II collagen by cathepsin K in equine articular cartilage. Osteoarthr. Cartil. 25, 2119–2126 (2017).

- Almeida, P. C. et al. Cathepsin B activity regulation. Heparin-like glycosaminogylcans protect human cathepsin B from alkaline pH-induced inactivation. J. Biol. Chem. 276, 944–51 (2001).
- Mehraban, F., Tindal, M. H., Proffitt, M. M. & Moskowitz, R. W. Temporal pattern of cysteine endopeptidase (cathepsin B) expression in cartilage and synovium from rabbit knees with experimental osteoarthritis: gene expression in chondrocytes in response to interleukin-1 and matrix depletion. *Ann. Rheum. Dis.* 56, 108–15 (1997).
- Ben-Aderet, L. *et al.* Detecting cathepsin activity in human osteoarthritis via activity-based probes. *Arthritis Res. Ther.* 17, 69 (2015).
 Charni-Ben Tabassi, N. *et al.* The type II collagen fragments Helix-II and CTX-II reveal different enzymatic pathways of human cartilage collagen degradation. *Osteoarthr. Cartil.* 16, 1183–91 (2008).
- Montaser, M., Lalmanach, G. & Mach, L. CA-074, but not its methyl ester CA-074Me, is a selective inhibitor of cathepsin B within living cells. *Biol. Chem.* 383, 1305–8 (2002).
- 52. Drake, F. H. et al. Cathepsin K, but not cathepsins B, L, or S, is abundantly expressed in human osteoclasts. J. Biol. Chem. 271, 12511-6 (1996).
- 53. Goto, T., Yamaza, T. & Tanaka, T. Cathepsins in the osteoclast. J. Electron Microsc. (Tokyo). 52, 551-8 (2003).
- Gelb, B. D., Shi, G. P., Chapman, H. A. & Desnick, R. J. Pycnodysostosis, a lysosomal disease caused by cathepsin K deficiency. Science 273, 1236–8 (1996).
- Leung, P., Pickarski, M., Zhuo, Y., Masarachia, P. J. & Duong, L. T. The effects of the cathepsin K inhibitor odanacatib on osteoclastic bone resorption and vesicular trafficking. Bone 49, 623–35 (2011).
- Zhuo, Y., Gauthier, J.-Y., Black, W. C., Percival, M. D. & Duong, L. T. Inhibition of bone resorption by the cathepsin K inhibitor odanacatib is fully reversible. *Bone* 67, 269–80 (2014).
- 57. Neutzsky-Wulff, A. V. et al. Alterations in osteoclast function and phenotype induced by different inhibitors of bone resorption implications for osteoclast quality. BMC Musculoskelet. Disord. 11, 109 (2010).
- Brage, M., Abrahamson, M., Lindström, V., Grubb, A. & Lerner, U. H. Different cysteine proteinases involved in bone resorption and osteoclast formation. *Calcif. Tissue Int.* 76, 439–47 (2005).
- 59. Everts, V. et al. Osteoclastic bone degradation and the role of different cysteine proteinases and matrix metalloproteinases: differences between calvaria and long bone. J. Bone Miner. Res. 21, 1399–408 (2006).
- 60. Henriksen, K. *et al.* Degradation of the organic phase of bone by osteoclasts: a secondary role for lysosomal acidification. *J. Bone Miner. Res.* **21**, 58–66 (2006).
- 61. Catterall, J., Dewitt Parr, S., Fagerlund, K. & Caterson, B. CTX-II is a marker of cartilage degradation but not of bone turnover Abstract. Osteoarthr. Cartil. 21, S77 (2013).
- 62. Bay-Jensen, A.-C. *et al.* The response to oestrogen deprivation of the cartilage collagen degradation marker, CTX-II, is unique compared with other markers of collagen turnover. *Arthritis Res. Ther.* **11**, R9 (2009).
- Mackie, E. J., Ahmed, Y. A., Tatarczuch, L., Chen, K.-S. & Mirams, M. Endochondral ossification: how cartilage is converted into bone in the developing skeleton. *Int. J. Biochem. Cell Biol.* 40, 46–62 (2008).
- Karsdal, M. A., Genovese, F., Madsen, E. A., Manon-Jensen, T. & Schuppan, D. Collagen and tissue turnover as a function of age: Implications for fibrosis. J. Hepatol. 64, 103–9 (2016).
- Mort, J. S. et al. Early cathepsin K degradation of type II collagen in vitro and in vivo in articular cartilage. Osteoarthr. Cartil. 24, 1461–9 (2016).
- 66. Combet, C., Blanchet, C., Geourjon, C. & Deléage, G. NPS@: network protein sequence analysis. Trends Biochem. Sci. 25, 147-50 (2000).
- 67. Gefter, M. L., Margulies, D. H. & Scharff, M. D. A simple method for polyethylene glycol-promoted hybridization of mouse myeloma cells. *Somatic Cell Genet.* **3**, 231–6 (1977).
- 68. Hui, W., Rowan, A. D., Richards, C. D. & Cawston, T. E. Oncostatin M in combination with tumor necrosis factor alpha induces cartilage damage and matrix metalloproteinase expression *in vitro* and *in vivo*. *Arthritis Rheum*. **48**, 3404–18 (2003).
- 69. Hao, J. L. et al. Effect of galardin on collagen degradation by Pseudomonas aeruginosa. Exp. Eye Res. 69, 595-601 (1999).
- Bossard, M. J. et al. Proteolytic activity of human osteoclast cathepsin K. Expression, purification, activation, and substrate identification. J. Biol. Chem. 271, 12517–24 (1996).
- Sørensen, M. G. et al. Characterization of osteoclasts derived from CD14+ monocytes isolated from peripheral blood. J. Bone Miner. Metab. 25, 36–45 (2007).
- 72. Henriksen, K., Karsdal, M. A., Taylor, A., Tosh, D. & Coxon, F. P. Generation of human osteoclasts from peripheral blood. *Methods Mol. Biol.* 816, 159–75 (2012).
- 73. Schindelin, J. et al. Fiji: an open-source platform for biological-image analysis. Nat. Methods 9, 676-82 (2012).

Acknowledgements

We acknowledge the funding by the Danish Science Foundation ("Den Danske Forskningsfond"). H.L. is supported by Marie Curie Initial Training Networks (Euroclast, FP7-People-2013-ITN: #607446). A.K. is supported by Marie Curie Initial Training Networks (BonePain, H2020-MSCA-ITN-2014: #642720).

Author Contributions

H.L., K.H. and C.S.T. conceived of and designed the study. H.L. developed and validated the GPDPLQ₁₂₃₇ ELISA, carried out biomarker measurements, and analysed the data. A.K. carried out BEX experiments and rat tissue isolation. A.D. performed immunohistochemical stainings. M.H.D. provided the human peripheral blood for osteoclast cultures. All authors contributed to data interpretation. H.L. wrote the manuscript. All authors critically reviewed the manuscript for intellectual content. All authors read and approved the final manuscript.

Additional Information

Supplementary information accompanies this paper at https://doi.org/10.1038/s41598-019-39803-0.

Competing Interests: H.L. performs research at Nordic Bioscience. A.K., M.A.K., Y.L., Y.H., T.M.J., A.C.B.J., C.S.T. and K.H. are employees of Nordic Bioscience. M.A.K., T.M.J., A.C.B.J. and K.H. hold stocks in Nordic Bioscience. Nordic Bioscience is a company involved in the discovery and development of biochemical biomarkers. A.D. and M.H.D. declare that they have no competing interests. The funders had no role in study design, data collection and analysis, decision to publish, or preparation of the manuscript.

Publisher's note: Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.

Open Access This article is licensed under a Creative Commons Attribution 4.0 International License, which permits use, sharing, adaptation, distribution and reproduction in any medium or format, as long as you give appropriate credit to the original author(s) and the source, provide a link to the Creative Commons license, and indicate if changes were made. The images or other third party material in this article are included in the article's Creative Commons license, unless indicated otherwise in a credit line to the material. If material is not included in the article's Creative Commons license and your intended use is not permitted by statutory regulation or exceeds the permitted use, you will need to obtain permission directly from the copyright holder. To view a copy of this license, visit http://creativecommons.org/licenses/by/4.0/.

© The Author(s) 2019